

Setup & Start Up Sequence

1. Remove the unit from the transport case (maintain vertical orientation).




2. Place the unit on a flat surface capable of supporting 42lbs (19kg).



3. Connect the main power cable to the unit and into an available electrical outlet (110V or 240V).



4. Apply power to the unit by switching the Main Power Switch on the back of the unit.

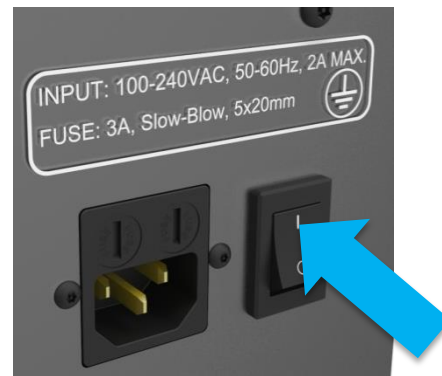

5. Power up the unit by pressing the Power Button on the front of the unit.



6. The ChipCHECK unit and software will initiate automatically. The unit is ready for use once the blue light is solid and the Welcome Screen is displayed.

